




PRODUCT INFORMATION

The Aeroneb® Go nebuliser is a fast, efficient, simple-to-use nebuliser developed for the millions of patients worldwide who require respiratory therapy in and away from home.

The Aeroneb Go can be used by paediatric as well as adult patients and is ideal for conditions such as asthma and chronic obstructive pulmonary disease (COPD). The Aeroneb Go aerosolises clinician-prescribed solutions for inhalation that are approved for use with a general purpose nebuliser.



Designed with convenience in mind and to help users maintain normal daily schedules, it is ultra-lightweight, portable, silent and is hassle free to set up and operate.

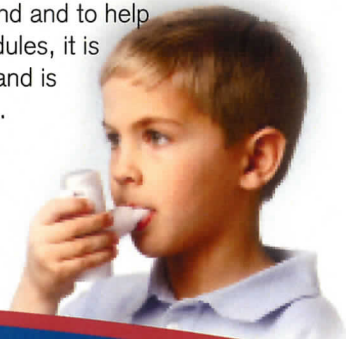

BENEFITS OF THE AERONEB GO


- **Fast:** Rapid treatment times (approx 7.5 minutes for a 3ml dose)
- **Lightweight:** 60g handset, fits comfortably into hand
- **Discreet:** Virtually silent
- **Effortless assembly and disassembly:** Ease of handling
- **Simple to clean:** Five straight forward methods recommended
- **Efficient:** Delivers a high-quality respirable aerosol. Aerosolises down to the last drop, minimises medication waste.
- **Ease of Use:** Simply fill the medication cup and press the On/Off button on controller
- **Convenient:** Operates off AA batteries or a mains electricity supply
- **Trouble free filling of medication cup:** Large opening (about the size of a €2 coin) with a simple twist off cap
- **Portable:** comes with a small carry case.

The Aeroneb Go incorporates Aerogen'd OnQ™ technology used in hospitals worldwide. OnQ™ technology is based on an electronic micropump that delivers aerosolised drugs to the lungs. It uses no compressed air, generates negligible heat and works with drugs in solution or suspension.


TWO POWER OPTIONS



1. AC/DC Powered

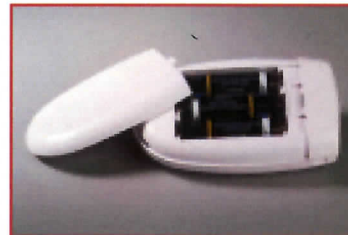

2. 3 AA Batteries

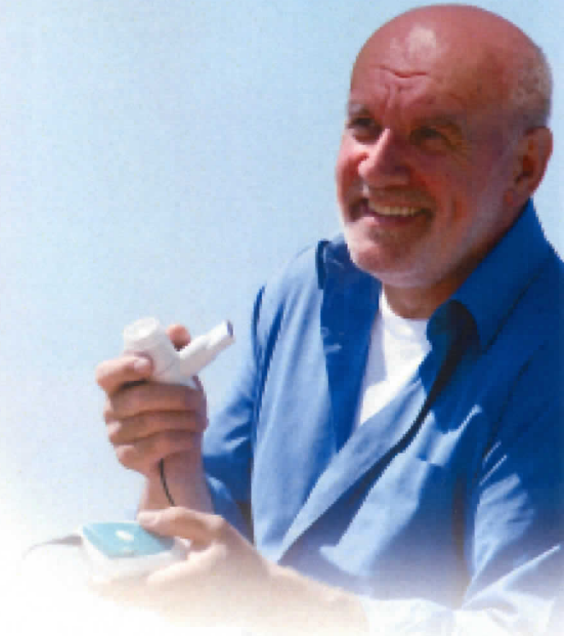

CLEANING OPTIONS

- Boil for up to 20 minutes
- Rinse under hot water
- Clean with a solution of warm water and a mild detergent
- Disinfect in a vinegar solution
- Sterilise in a baby bottle steam steriliser



SPECIFICATIONS

Weight	Nebuliser Unit 60g
Sound level	<35dBA at 1 metre
Medication cup capacity	6ml
Residual volume typically	Typically <0.02ml
Average flow rate	0.4ml/min
Average particle size	3.6µm MMAD

LIFESTYLE NEBULISER

- Fast
- Simple
- Lightweight
- Efficient

For more product information, including accessories and parts, contact your local supplier:

LOCAL DISTRIBUTOR

DISTRIBUTED IN AUSTRALIA BY:

**Tag
Medical**

www.tagmedical.com.au
P.O. Box 428, Canterbury, Vic. 3126
T: (03) 9888 5111 F: (03) 9836 8840

Aerogen Ireland Ltd,
Galway Business Park, Dangan, Galway, Ireland
Phone: +353 91 502550
Fax: +353 91 584639

www.aerogen.com

CE
0050

Aerogen and Aeronneb are registered trademarks of Aerogen, Ireland Ltd. PM/022


Aerogen®


Aeroneb® Go
Micropump Nebuliser




Aerogen®